


USED ON


IF IN DOUBT ASK

DO NOT SCALE

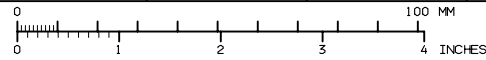


PART MUST BE FREE FROM BURRS

DRG. NO. NPS-0047-I

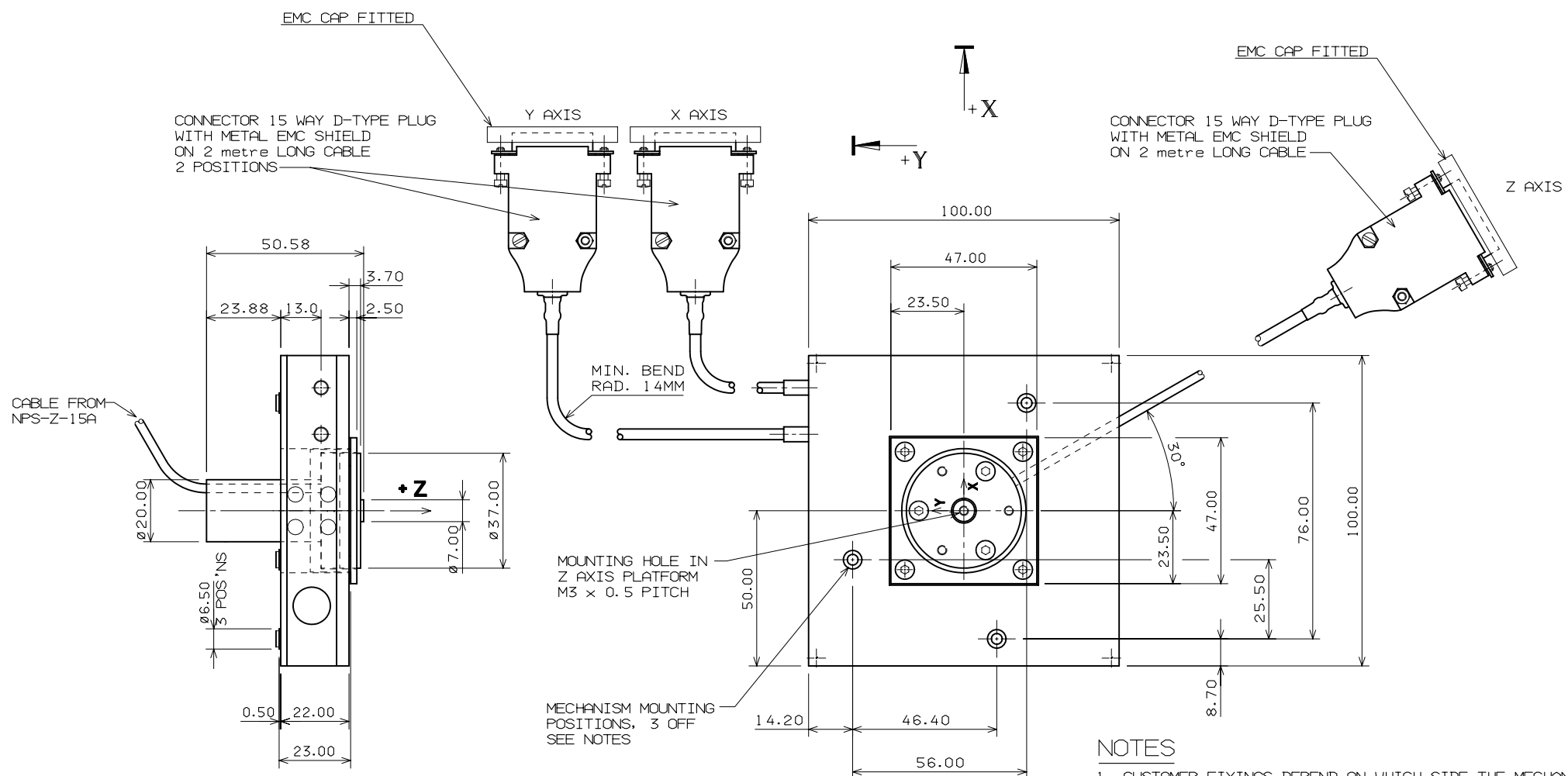

NOTES

1. CUSTOMER FIXINGS DEPEND ON WHICH SIDE THE MECHANISM IS TO BE MOUNTED FROM.
2. TO MOUNT MECHANISM FROM ABOVE, 3 OFF M3 CAP HEAD SCREWS ARE REQUIRED.
3. TO MOUNT MECHANISM FROM UNDERNEATH, 3 OFF M4 SCREWS ARE REQUIRED.
4. ALL FIXING POSITIONS REQUIRE A COUNTERBORE OF Ø5.25 x 1.5mm DEEP TO LOCATE MECHANISM.

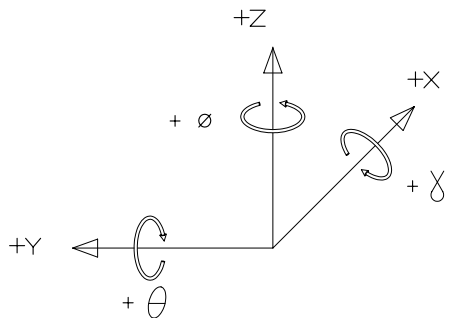

ISSUE	ECR No	DESN	DRAWN	CHKD	REA	APP'D	DATE
1		YX/JR	GMJ				

ALL DIMENSIONS IN MM	
TOLERANCES UNLESS OTHERWISE STATED	
1 DECIMAL PLACE	±0.25
2 DECIMAL PLACES	±0.15
TRUE POSITIONAL ZONE DIA.	0.1
ANGULAR (DEGREES)	±0.5

© COPYRIGHT 1999 QUEENSGATE INST. LTD.	QUEENSGATE INSTRUMENTS LIMITED WILLOUGHBY ROAD, BRACKNELL, BERKSHIRE, RG12 8FB.	
MATERIAL	TITLE INSTALLATION NPS-XY-100/NPS-Z-15A	
FINISH	SCALE 1-1	DRG. SIZE A2
	DRAWING NUMBER NPS-0047-I	ISSUE 1
		SHEET OF 1

